

SMACC

Sleep Manipulation And Communication Clickything — manual

SMACC

None

Table of contents


1. SMACC	4
1.1 What SMACC does	4
1.2 Get started	5
1.3 Used by	5
2. Getting started	7
2.1 Installation	7
2.2 SMACC files	9
3. Running a session	11
3.1 Overview	11
3.2 Audio cues	0
3.3 Visual cues	0
3.4 Biocals	0
3.5 Dream reports & surveys	0
3.6 Intercom & chat	0
3.7 Markers & port codes	0
4. Devices & routing	0
4.1 Audio routing	0
4.2 Compatible devices	0
5. Volume & latency	0
5.1 Volume	0
5.2 Latency	0
5.3 What contributes	0
5.4 The marker marks the software event, not the photons	0
5.5 Measured: a typical Windows rig	0
5.6 The Low / High setting	0
5.7 Visual cues: BlinkStick vs Hue	0
5.8 Measure it on your own rig	0
5.9 If you want it lower	0
6. EEG Annotator	0
6.1 Supported recordings	0
6.2 Viewing	0
6.3 Annotating	0
6.4 Sleep staging	0
6.5 Session log overlay	0
6.6 Multiple raters	0

6.7	Blind-rater mode	0
6.8	For developers	0
7.	Troubleshooting	0
7.1	Common problems	0
7.2	If SMACC crashes	0
7.3	Reporting a problem	0
8.	Release notes	0
8.1	0.1.0	0
9.	Reference	0
9.1	File formats	0
9.2	SMACC file (<code>.smacc</code>)	0
9.3	Preferences file (<code>preferences.yaml</code>)	0
9.4	Session log (<code>.log</code>)	0
9.5	BIDS export (<code>events.tsv</code> + <code>sidecar</code>)	0
9.6	Annotations sidecar	0
10.	Contributing	0
10.1	Conventions	0
10.2	Commit and pull-request style	0
10.3	Development	0
10.4	Building the executable	0
10.5	Building the docs	0
10.6	Versioning	0
10.7	Project notes	0

1. SMACC

Sleep Manipulation And Communication Clickything (SMACC) is a Windows desktop app for running sleep and dream studies — presenting cues to a sleeping participant, communicating with them, marking events on the EEG, and collecting dream reports.

In a dream-engineering session the experimenter has to deliver a precisely-timed cue to a sleeping participant, mark each event so it lines up with the EEG, talk with the sleeper, and collect a dream report — all from a dark control room, late at night, on whatever hardware the lab has. SMACC is the clickable control surface for that work: one window per job, glanceable in the dark, with a volume safety cap and explicit device routing so a misclick can't blast or mislead a sleeping participant.

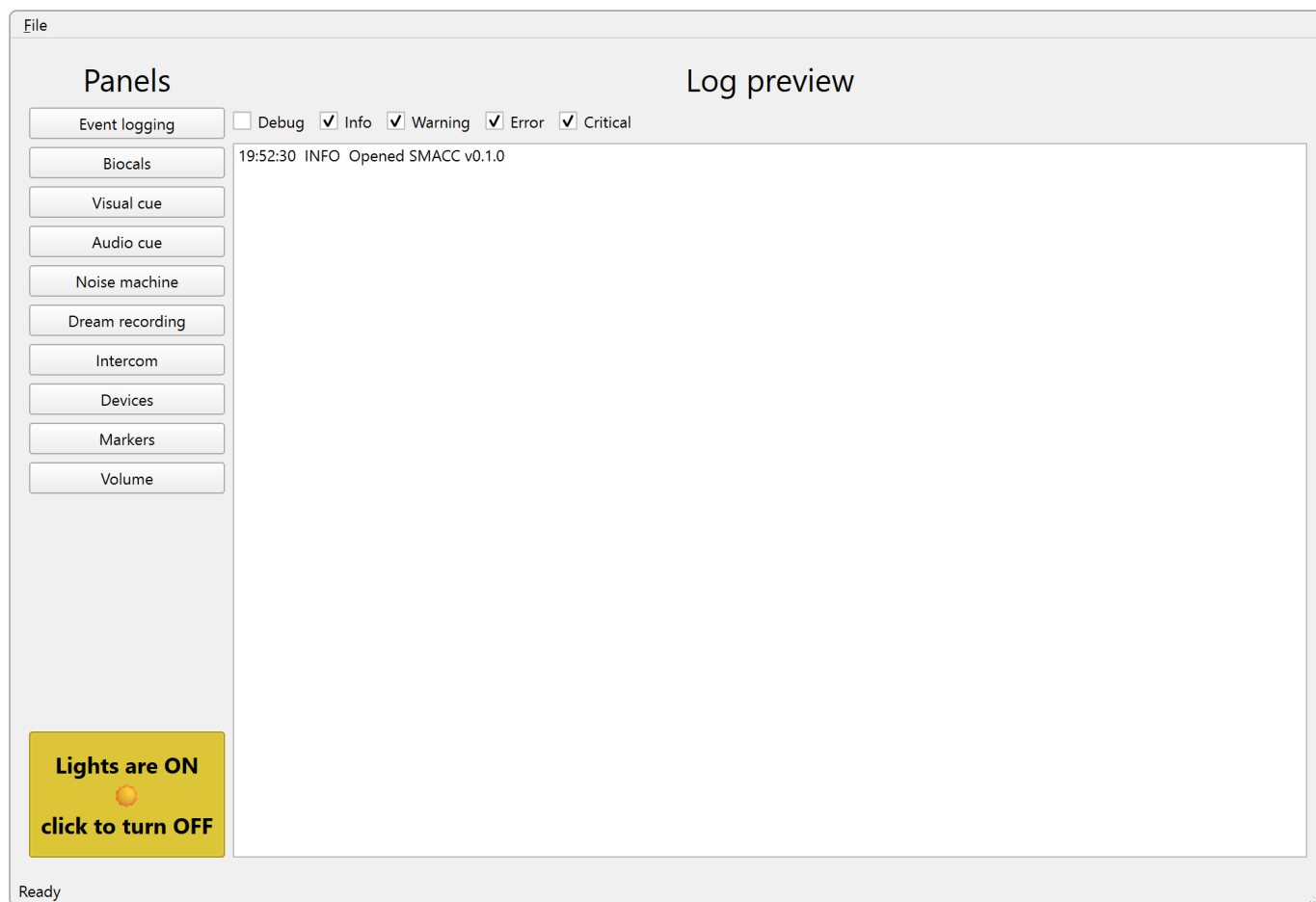
[Download SMACC](#)
[Installation guide](#)

The button downloads the installer (`SMACC-Setup.exe`) for the latest version, which installs per-user (no admin rights) and runs on 64-bit Windows 10 or later. Prefer a single file with no install? A portable `SMACC.exe` ships with every release — see [Installation](#), which also covers the Windows SmartScreen warning on the unsigned download and how to get older versions.

Working offline? [Download these docs as a single PDF](#).

1.1 What SMACC does

- Trigger **audio cues** (and design them in-app with the Audio Cue Designer)
- Trigger **visual cues** on a BlinkStick or Philips Hue light
- Play masking **background noise**

- Run **biocalibrations** as timed, marked tasks
- Record **dream reports** and administer **surveys**
- **Talk, listen, and type** with the participant over an intercom
- Mark events with **EEG port codes** over LSL or a hardware TTL trigger
- Review and score recordings in the **EEG Annotator**
- Save a detailed **event log** for every run

1.2 Get started

Getting started

- [Installation](#) — download and run SMACC.
- [SMACC files](#) — create and reuse a study's configuration (`.smacc`).

Running a session

- [Overview](#) — the Launcher, the Session window, and the tools.
- [Audio cues](#) · [Visual cues](#) · [Biocals](#) · [Dream reports & surveys](#) · [Intercom & chat](#) · [Markers & port codes](#)

Devices, volume & timing

- [Audio routing](#) · [Compatible devices](#) · [Volume & latency](#)

After the night

- [EEG Annotator](#) — review and score recorded EEG.
- [Troubleshooting](#) · [Reference](#) · [Contributing](#)

1.3 Used by

SMACC is used in dream engineering research, including by:

- [Ken Paller's Cognitive Neuroscience Lab](#) at Northwestern University
- Torres-Platas et al. (2026). Intentional lucid dreaming with a transformative learning agenda. *Research Square* doi:[10.21203/rs.3.rs-8745420/v1](#)
- Konkoly et al. (2026). Using real-time reporting to investigate visual experiences in dreams. *J Cogn Neurosci* doi:[10.1162/jocn.a.107](#)
- Morris et al. (2026). Inducing lucid dreaming based on a contemplative practice of compassion. *Brain Sci* doi:[10.3390/brainsci16030315](#)
- Konkoly et al. (2026). Creative problem-solving after experimentally provoking dreams of unsolved puzzles during REM sleep. *Neurosci Conscious* doi:[10.1093/nc/niaf067](#)
- Konkoly et al. (2025). Investigating dreams by strategically presenting sounds during REM sleep to reactivate waking experiences. *Neuropsychologia* doi:[10.1016/j.neuropsychologia.2025.109229](#)
- Morris et al. (2025). Lucid dreaming of a prior virtual-reality experience with ego-transcendent qualities: A proof-of-concept study. *Neurosci Conscious* doi:[10.1093/nc/niaf017](#)
- Mundt et al. (2024). Treating narcolepsy-related nightmares with cognitive behavioural therapy and targeted lucidity reactivation: A pilot study. *J Sleep Res* doi:[10.1111/jsr.14384](#)
- Wolk et al. (2024). Lucid dreams from reactivating mindfulness during REM sleep: A pilot study. *Int J Dream Res* doi:[10.11588/ijodr.2024.2.98233](#)
- [Michelle Carr's Dream Engineering Lab](#) at the University of Montreal and the Center for Advanced Research in Sleep Medicine
- Jafarzadeh Esfahani et al. (2024). Highly effective verified lucid dream induction using combined cognitive-sensory training and wearable EEG: A multi-centre study. *bioRxiv* doi:[10.1101/2024.06.21.600133](#)

SMACC is free software released under the [GPL-3.0-or-later](#) license.

2. Getting started

2.1 Installation

Download SMACC

The button above always downloads the installer for the **latest** version — double-click `SMACC-Setup.exe` and click through. The installer:

- installs SMACC **per-user** — no administrator rights or IT involvement needed;
- adds a **Start menu** entry (and, if you opt in, a desktop shortcut);
- associates `.smacc` files, so double-clicking one opens the **Start a session** dialog with it preselected;
- installs the **EEG Annotator** — the post-hoc [EEG viewer/annotator](#). Included by default; uncheck it during setup to skip the heavyweight EEG component (which bundles the MNE library), or re-run the installer later to add or remove it;
- registers an uninstaller — remove SMACC from **Settings > Apps** like any other program. Uninstalling never touches your data: SMACC files, recordings, and logs under `~/SMACC` (or your data directories) all stay put. The uninstaller reminds you of this (and of the folder's location) when it finishes, so finding the folder afterwards is expected — delete it manually if you no longer need the data.

Installing a newer version over an existing one upgrades it in place, keeping whichever components you had chosen.

To get an older or specific version, browse the [releases page](#): pick the release you want, open its *Assets* dropdown, and download from there. (The version switcher on this documentation site switches the *docs* only — older copies of SMACC itself come from the releases page.)

System requirements

SMACC runs on 64-bit Windows 10 or later.

Windows SmartScreen — “Windows protected your PC”

SMACC isn't code-signed yet, so when you run the installer Windows SmartScreen may show a blue “**Windows protected your PC**” box. This is expected for any new, unsigned program. To proceed anyway, click **More info**, then **Run anyway**.

Administrator privileges and the UAC prompt

For some features you will need to open SMACC with Administrator privileges (right-click the Start menu entry and select **Run as administrator**). Windows then shows a **User Account Control (UAC)** prompt asking whether to allow the app to make changes — click **Yes** to continue. For everyday use (audio cues, dream reports, LSL markers) you can run SMACC normally, without administrator rights.

2.1.1 Portable SMACC.exe (no install)

Every release also ships the same app as a single portable `SMACC.exe` — download it, double-click it, no installation. This suits USB-stick deployment, machines where nothing may be installed, and quick tests. The portable build skips the installer's conveniences (Start menu entry, uninstaller, the `.smacc` file association) — you can associate `.smacc` files from the Launcher's File menu instead (see below). The **EEG Annotator** is likewise available portable as `SMACC-EEG.exe` — put it **next to** `SMACC.exe` and the Launcher's *EEG Annotator* button lights up (it also runs standalone).

For IT departments

The default install is per-user (under `%LOCALAPPDATA%\Programs\SMACC`, no elevation). On managed machines where per-user installs are blocked — e.g. AppLocker policies on `%LOCALAPPDATA%` — run the same installer machine-wide instead: `SMACC-Setup.exe /ALLUSERS` (elevates, installs to Program Files).

2.1.2 Updating SMACC

From the Launcher, **File** > **Check for updates...** asks GitHub whether a newer release exists and, if one does, offers to open the download page in your browser — run the new installer and it upgrades the existing install in place. SMACC never checks on its own: lab machines are often offline, and studies usually pin one version for their whole run, so checking is always an explicit click. If you use the portable `SMACC.exe`, download the new one and replace your old copy.

2.1.3 After installing

By default SMACC stores everything under `~/SMACC`. To use a different location, set the environment variable `SMACC_DIRECTORY` to the directory you want; SMACC creates it and its subfolders on first run.

The rest of setup is covered on the relevant pages:


- **Your study configuration** lives in a portable [SMACC file](#). SMACC seeds a `default.smacc` and opens it when you don't pick another, so it works out of the box. The installer associates `.smacc` files; with the portable build, enable or repair the association from **File** > **Associate .smacc files (Windows)** in the Launcher.
- **Audio cues** go in the data directory's `cues/` folder, alongside the seeded `demo-*` cues — see [Audio cues](#).
- **Dream-report surveys** are managed from the Dream recording panel — see [Dream reports & surveys](#).
- **A recording mic**: to record dreams, bind your mic to **Bedroom mic 1** in the **Devices** window (the **Panels** column) — see [Audio routing](#).

2.2 SMACC files

A **SMACC file** (`.smacc`) is the study's configuration: cue files, volumes, noise, visual cues, biocals, survey presets, event markers, the display choices that apply to a session (always-on-top and log-preview levels), and the **data directory** where runs are written. It plays the same role for a SMACC study that a montage file plays for an EEG setup — configure it once, then reuse it session after session.

A SMACC file is *not* a session file. It is not tied to any particular run: opening one **starts a new session** with that configuration, and you can start as many sessions from the same file as you like. Sessions never rewrite the SMACC file they were started from — what happened in a run is recorded in that run's [session log](#), not in the SMACC file.

It is a single, portable file (plain YAML you can read and edit). Cue/noise files and the data directory are stored **relative** to the SMACC file when they sit beside it — so a self-contained study folder stays valid if you move, copy, or zip it — and **absolute** when they point elsewhere. For the exact on-disk format, see the [SMACC file reference](#).

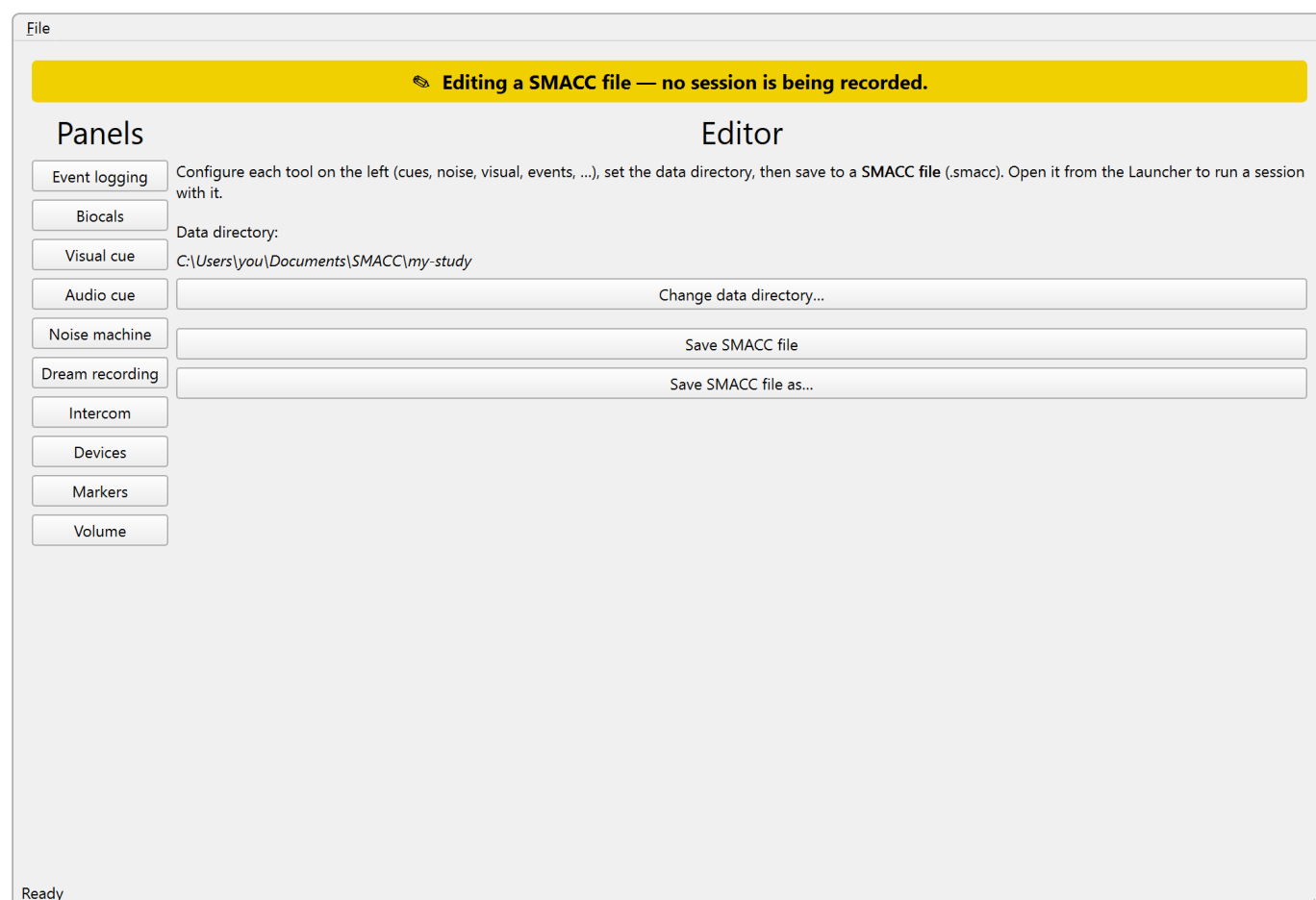

2.2.1 Creating a SMACC file

There are two ways to make one:

- **From scratch, in the Editor.** In the [Launcher](#), click **Editor...**. In the dialog, choose **New SMACC file** for a blank configuration (or pick an existing file to edit). Configure each panel on the left, set the data directory, and **Save SMACC file as...** to a new `.smacc` anywhere you like. The Editor never records a run — it only authors configuration.
- **From a live session.** In a running **Session**, use **File > Save SMACC file as...** to snapshot the session's *current* settings — including any mid-session tweaks (volumes, device routing, event-code edits) — to a new SMACC file. Handy when you've tuned things live and want tonight's setup to be tomorrow's starting point.

SMACC ships a `default.smacc` in your SMACC directory as a working example. It is treated as a read-only template — saving over it redirects to Save-As — so there is always a known-good configuration to fall back on.

You can keep SMACC files anywhere: for instance one per participant (`peter.smacc` , `paul.smacc` , ...) or one per protocol, each pointing at whatever data directory you like (they can share one).

2.2.2 The data directory

Every SMACC file names a **data directory** — where its runs are written. Each session started from the file gets its own timestamped folder under that directory (e.g. `smacc-20260607-223015/`) holding the run's `.log` , dream-report recordings, survey responses, and any exports. Set it in the Editor (**Change data directory...**); choose it deliberately, since it determines where a night's data ends up.

2.2.3 Opening a SMACC file

Opening a SMACC file starts a new session configured by it:

- **Double-click the file** (Windows build): SMACC opens the **Start a session** dialog with that file preselected, so you confirm the subject/session and click **OK** to start. Enable or repair the file association any time from **File > Associate .smacc files (Windows)** in the Launcher.
- **From the Launcher**: click **Session...**, then pick the file in the **Start a session** dialog (recent files are listed; **Browse...** finds any other) and click **OK**.
- **From a terminal**: `SMACC path/to/file.smacc`.

3. Running a session

3.1 Overview

SMACC opens to its **Launcher** — a list of the SMACC tools. Click a tool and SMACC asks for the **SMACC file** to use, then opens it. To run a night you open a live **Session**. This page orients you; each tool then has its own page.

The Launcher's tools (its title bar reads **SMACC**):

- **Session...** — run a live **Session** for collecting data, the interface for running a night. SMACC prompts for a [SMACC file](#) first.
- **Editor...** — create or edit a [SMACC file](#) in the **Editor**, without recording anything.
- **Analyzer** — inspect a past session.
- **Audio Cue Designer** — build a tone cue and export it as a WAV.
- **EEG Annotator** — review and annotate a recorded EEG (see [EEG Annotator](#)); disabled when its optional component is not installed.

The **Session** and the **Editor** are the same window in two modes — running a night versus authoring a SMACC file.

File



SMACC

Session...

Editor...

Analyzer

Audio Cue Designer

EEG Annotator

v0.1.0 — [documentation](#)